



← Retro-fit any existing X-Cel X-Ray unit

The PDRW is cassette sized, giving you the ability to retro-fit any X-Cel unit. 20/20 Imaging also offers new x-ray units for practices without existing x-ray, or those who are simply looking to upgrade their older x-ray equipment.

All X-Cel units including: MB700 (left), P75/P700, FB700, BD715, Wall-Mount, Mobile

↓ New X-Ray

2020's PXS710-D (right) offers a High Frequency digital x-ray solution with full software integration including anatomical programming. The PXS710 uses ½ the radiation dose, providing twice the efficiency, and enhanced, sharper image quality over standard frequency with a 0.5mm focal spot. Paired with the new PDRW, it offers imaging excellence to the podiatry office.



P-DRW Wireless Panel



Taking digital radiography to the next level with our Wireless Flat panel technology is 2020 Imaging's cassette-sized PDRW. Providing exceptional image quality, paired with advanced software features and functionality. The PDRW panel uses enhanced Gadolinium technology and has many unique features such as device-drop notification emails.

The PDRW is ideal for use with both new x-ray units as well as a retrofit solution. The versatile panel utilizes the Opal-RAD Acquire software for image capture, manipulation, and sending. 2020's acquire software offers useful tools to enhance the image acquisition process, thus keeping your office in the forefront of x-ray technology.

SPECIFICATIONS

- Warranty: 3 years
- Wireless connectivity
- Cassette sized
- Battery Life: > 3.5 hrs
- 7.93 lbs
- Pixel Matrix: 2304x2800
- Off-center Imaging
- Dimensions: 15.12"x18.11"x0.6"
- A/D conversion: 14-bit



Customized Practice-Specific Software*

20/20 Imaging integrates the specialized Opal Podiatry tool-set for the Podiatric discipline offering advanced measurement and alignment tools. Every New or Retro-fit 2020 system comes complete with a set of advanced PACS tools and functionalities for your office.

*Available for an additional cost. Contact sales for details.




X-CEL X-RAY RETRO-FIT KITS

FULLY DIGITAL WIRELESS FLATPANEL SOLUTION



FB700

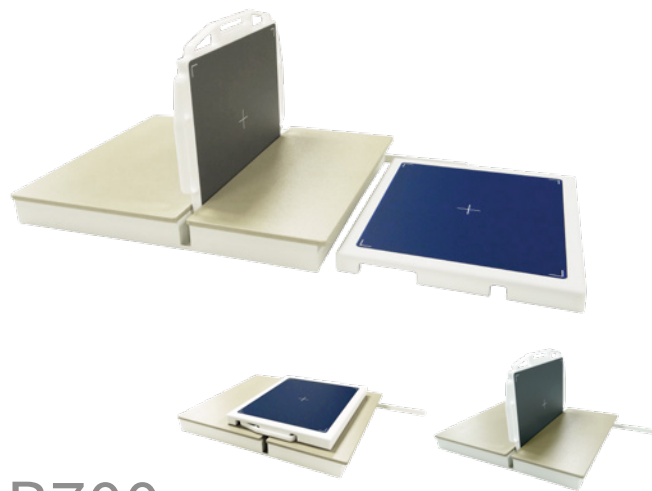

LB700

MB700

P700/P75

