

MOMENTUM-EQ

ALL-IN-ONE EQUINE WIRELESS DR SOLUTION

Introducing the new Momentum-EQ, the most lightweight, portable, all-in-one Equine DR solution available in it's class.



All-in-one Charging
Charge everything with one power cable


- Lightweight - under 18 lbs
Custom holders and design available on MO-EQ
- Built-in positioner



- Ergonomic custom panel and tablet holders
Accommodates all size hands for easy handling
- Rugged
Lightweight aluminum durable design
- Scratch resistant finish

10"x12" Premium Wireless DR Panel

Manufactured by Konica Minolta

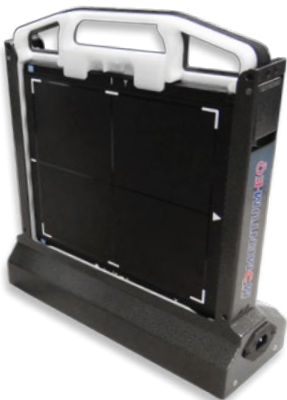

DR Panel Specifications

- Lightest panel & holder at 4lbs
- Scintillator: CsI (Cesium)
- Cassette sized
- Imaging area: 10" x 12"
- 16-bit AD Conversion
- Battery life: > 150 images/4hrs
- Rapid charge enabled
- Warranty: 3 years
(7 year capacitor)

- Custom tablet enclosure
Freestanding or hung on stall to eliminate the need to bend down or kneel
- All-in-one Software Solution
DR Console, Clinical Review, and Image Storage
- User-friendly interface
- Automated image processing
Produces an optimized image - every time

